
## PRODUCT DESCRIPTION

REGLOW LED FACE MASK

USER MANUAL

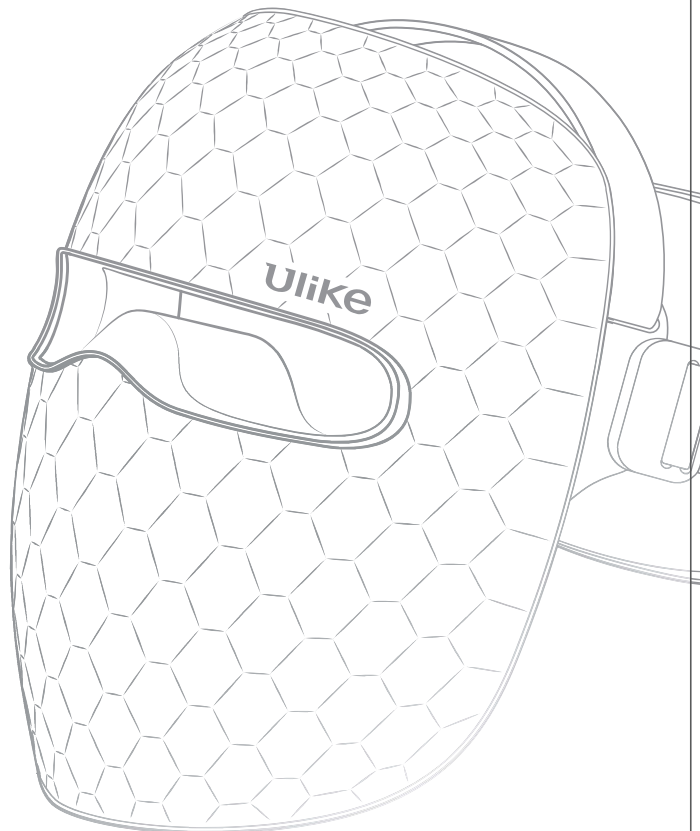Please read this user manual carefully before using the product.

ENGLISH ————— 02

Hebrew ————— 16

Thank you for choosing our product.  
To ensure safe and proper use, please  
carefully read this manual and follow  
its instructions. The images in this  
manual are for illustrative purposes  
only, which may differ from the actual  
product.

For any issues during use, please contact our customer service at:  
[seasupport@ulike.com](mailto:seasupport@ulike.com).

## I . Safety Precautions

### ● Important notes

- Please assemble and use the product in accordance with the User Manual.
- Please carefully read the User Manual and keep it for future reference.
- If you have any questions or concerns during use, please contact our online customer support.

❗ **Danger:** Improper use of this device may cause injury.

⚠ **Caution:** Improper use of this device may cause injury or property damage.

ⓘ **Warning:** Improper use of this device may lead to dangerous situations.

⚠ **Caution:** Improper use of this device may cause injury or property damage.

ⓘ **Warning**

Do not use the device under the following conditions. Failure to comply may cause accidents or harm to skin and body.

- If you have a pacemaker, defibrillator, or other implanted medical device.
- If you are under the influence of alcohol, drugs, or medication that may impair your judgment or ability to operate the device safely.
- If you have a reduced or impaired sensation of temperature or pain in the area of application.
- If you have an autoimmune disease.

※ **If you have had surgery within the past 2 weeks, consult your doctor before using the device. Before using this device, consult a medical professional if you have any of the following conditions. Using this device without proper consultation may cause harm to your skin and body.**

- An acute illness, infectious disease, malignant tumor, high blood pressure, blood disease, fever, or febrile illness.
- Sensory impairment due to severe peripheral circulatory dysfunction (e.g., diabetes), nerve damage, skin disease, or atopic dermatitis, or if you are pregnant or breast feeding
- A hospital outpatient status and taking medication, a pre-existing skin condition, such as an allergic constitution or contact dermatitis, or a concern about your skin when using this device.

Follow these precautions to avoid possible damage to your skin and body.

- Do not use the device for an extended period of time or over the same area more times than recommended.

- Stop using the device immediately if you experience discomfort or adverse effects, such as an allergic reaction, skin disease, allergic dermatitis, or sensitive skin.

- Follow the user manual instructions carefully. Failure to do so may result in adverse effects.

This device contains a rechargeable battery, which may overheat, catch fire, or explode if not handled properly. To prevent these risks, follow these safety guidelines.


- Do not charge, use, or place the device near an open flame and heating devices or in any high-temperature environments.
- Keep out of reach of children. Do not use on children or allow children to use the device, as this may cause accidents or injury.
- Avoid use on those who cannot express their intentions, and do not allow disabled individuals to use the device unattended as this may cause accidents.
- Use only as directed in this manual.
- Do not use near flammable substances such as alcohol, thinner, gasoline, nail polish remover, or hair spray as this may cause fire, electric shock, or injury.

**⚠ CAUTION**

- Do not damage, process, bend by force, pull, twist, or bundle the cord. Do not place heavy objects on the cord or allow it to get pinched in between other objects as this may damage the cord and cause a fire or electric shock.
- Do not use a damaged Type-C charging cable or loose socket to avoid potential electric shock, short circuit, or fire. Consult the distributor or after-sales service if repair is needed.
- Keep the main unit and charger dry; do not wash them with water.
- Do not use near water sources, or while in the bath or shower to avoid potential electric shock, short circuit, or fire.
- Do not disassemble, repair, or modify this device in any way as this may cause fire, electric shock, or injury. Consult the distributor or after-sales service if repair is needed. Use an adapter that complies with national standards and holds CQC/CCC certification, with a rated input voltage of DC 5V and output power ranging from 5W to 15W.
- Do not use the device while charging.
- Do not place the device near fire, under direct sunlight, or in humid places such as the bathroom. Do not place near heating devices as this may cause overheating, fire, or even explosion. Immediately stop charging in case of abnormal situations such as smoke or an unusual noise or odor.
- Do not use if you have any unusual skin condition such as a rash, irritation, or itching. Avoid excessive or continuous use of the device.
- Stop using the device and unplug it from the socket immediately if you notice any abnormal odor or unusual overheating. This may indicate damage to the device or a risk of fire.
- Before charging, make sure the power plug and socket are free from moisture to avoid potential electric shock, short circuit, or fire.
- Be sure to turn off the device when not in use or when being cleaned to avoid potential electric shock, short circuit, or fire.
- Do not use the device if it has been damaged, dropped, or subjected to a strong impact. This may affect the performance of the device or cause injury.
- Do not wipe the main unit with alcohol, thinner, gasoline, nail polish remover, or any solvents which could damage the surface of the device.

- Place the device on a stable and flat surface when not in use. Avoid dropping or subjecting the device to impact as this may cause malfunction or injury.
- If you see any abnormal skin reaction during or after use, stop using the device and consult a doctor.
- Do not look directly at the LED lights to prevent possible eye discomfort.
- In the event of a power outage, turn off the device and unplug it from the socket immediately. This will prevent accidents or damage to the device when the power is restored.
- Pay attention to your skin's condition while using the device. If you feel skin fatigue or any unusual sensations, stop using it immediately. If the condition persists the next day, consult a doctor as soon as possible.
- The appliances are not to be used by :
  - Persons suffering from sunburn
  - Persons under medical care for diseases that involve photosensitivity;
  - Persons receiving photosensitising medications
  - Persons having a first-degree relative with a history of melanoma;
  - when the appliance is used in a bathroom, unplug it after use since the
- Proximity of water presents a hazard even when the appliance is switched off;
- The appliance should be cleaned to avoid the accumulation of grease and other residues.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical,
- Sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or
- Instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

## II. Product description and accessories

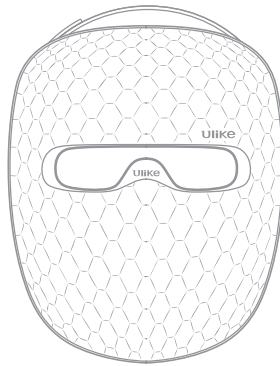
LED Mask


Rechargeable  
Remote


Type-C charging cable

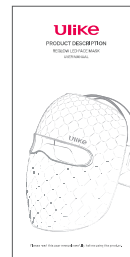User Manual
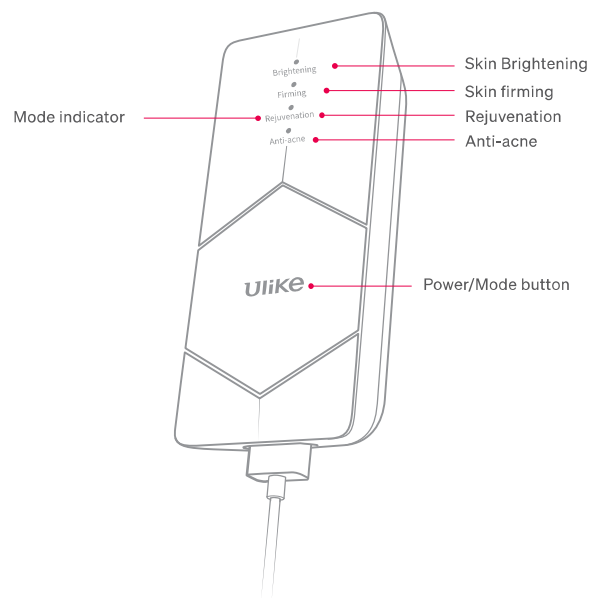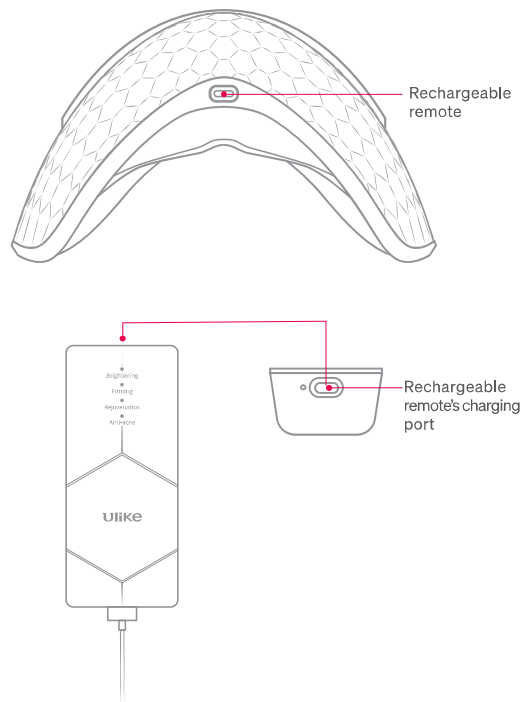
#### Charging the device

The device consists of a LED facemask and a rechargeable remote. You only need to charge the rechargeable remote, as it will power the main unit automatically.

To charge the rechargeable remote, connect the Type-C charging cable to the charging port on the remote. Plug the other end of the cable into a compatible power adapter (not included).

The charging indicator on the remote shows the battery status of the device:

- When the device is on, a flashing red light indicates that the battery is low and needs charging. While charging, the red light stays solid. When the battery is fully charged, the light turns green.

- When the device is off, you can press the power button to check the battery status. If the battery is low, the charging indicator will flash red eight times and then go out. This means you need to charge the device. If the battery is sufficient, the charging indicator will flash green once and then go out.



### III. How to use the device

#### Before using the device for the first time

##### Precautionary check


- Check the device and accessories. Take out the main unit and accessories from the box. Keep the packaging materials for future transportation, storage, and reference. Check the main unit and accessories for any signs of mechanical damage. If you find any damage, do not use the device, and contact the distributor or after-sales service as soon as possible.
- Perform a light sensitivity test.

##### Light sensitivity test

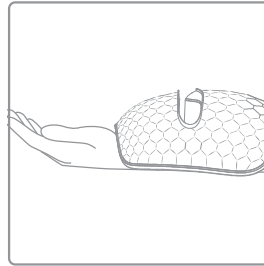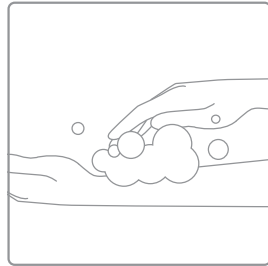
- Cleanse the inner arm with soap or a facial cleanser, pat it dry, and then expose it to the light from the LED mask for 3 minutes. Wait for 24 hours and check your skin for any unusual conditions, such as persistent redness, burning, or heat. If you do not have any adverse reaction, you can use the device on your face. However, if you have any rash, redness, itching, or other unusual reactions before, during, or after the test, stop using the device immediately.


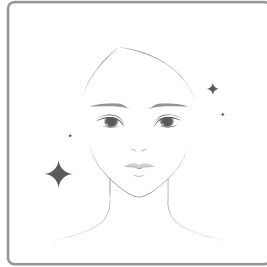
1. Remove your makeup and cleanse your face before use. Apply toner for basic moisturizing, but be cautious not to use skincare products that are light sensitive or contain whitening ingredients.

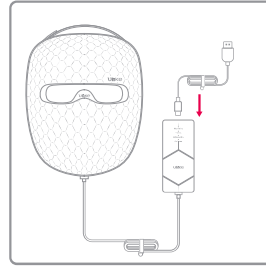2. Connect the mask and controller with the provided cable.

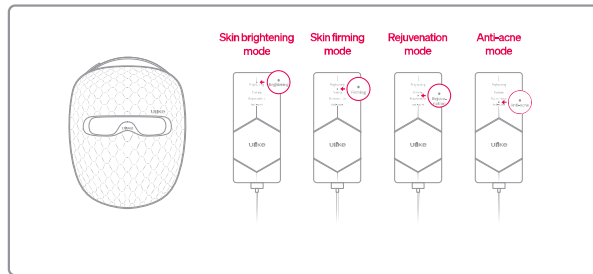3. This device features a memory function. When you turn it on, it will resume the mode used during the last session. Select a mode when using it for the first time. Press and hold the Power/Mode button on the controller for 1 to 2 seconds to turn on the device. The treatment will automatically start with Skin brightening mode. Press the Power/Mode button to cycle between Skin brightening, Skin firming, Rejuvenation, and Anti-acne modes.

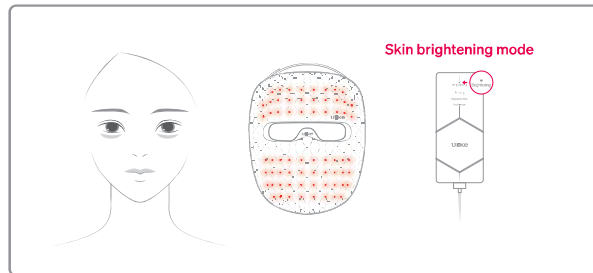
#### 4. Treatment

##### • Skin brightening mode

**Ideal for addressing the impact of late nights or for daily skin rejuvenation to counteract a sallow complexion.**

Tips: After makeup removal and cleansing, set the device to this mode. Following the light therapy, apply whitening skincare products for enhanced effects. Do not use light-sensitive or irritating products containing vitamin C, retinol (vitamin A), or similar substances before the light therapy.

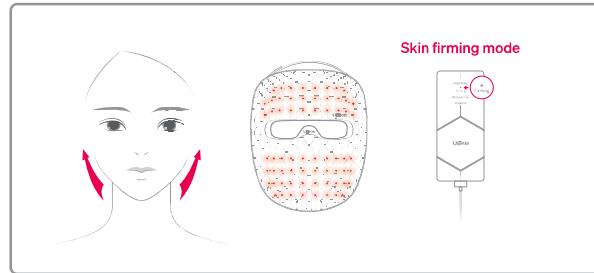
#### 4. Treatment

- Skin firming mode

**Ideal for fighting signs of aging or early aging.**

Tips: After makeup removal and cleansing, switch to this mode for treatment. After the light therapy, apply anti-aging products for boosted results.

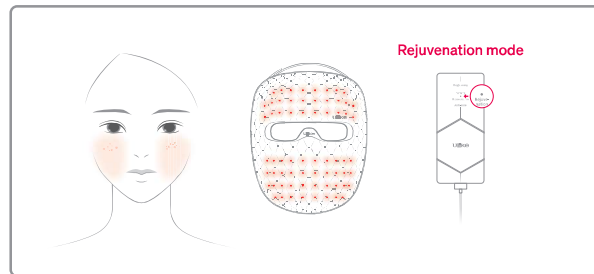
#### 4. Treatment

- Rejuvenation mode

**Ideal for sunburn recovery, seasonal sensitivity, allergy relief, or recovery after cosmetic surgery.**

Tips: After pampering your skin with a gentle, irritation-free cleanser, switch to this mode for treatment. During sensitive periods when the skin barrier is fragile, make sun protection a daily ritual.

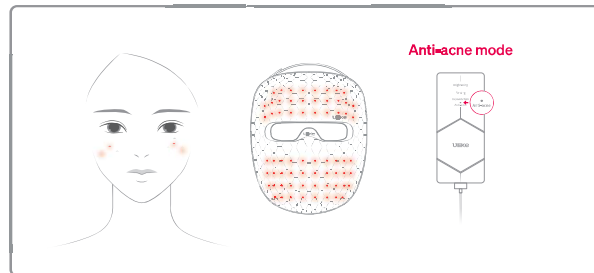
#### 4. Treatment

- Anti-acne mode

**Ideal for creating a perfect balance between moisture and oil, refining clogged pores, or reducing acne.**

Tips: After makeup removal and cleansing, apply a light moisturizer before activating this mode. If your skin feels dry after the session, boost hydration with extra moisturizing products. It is also advisable to maintain a healthy lifestyle for more noticeable results.

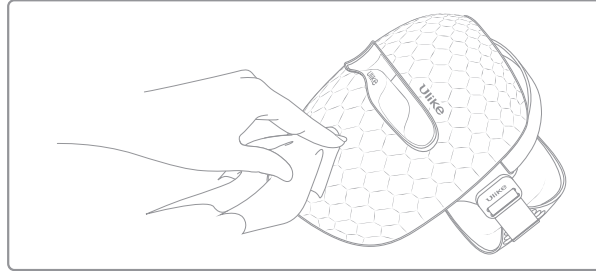**5. After treatment**

After treatment, place the device back into the packaging box or storage bag. If it gets dirty, clean it with a slightly damp soft cloth or sanitizing wipe.

**Do not use alcohol, thinner, gasoline, nail polish remover, or any solvents to clean the device**

**6. Frequency of use**

Use up to 6 times a week, choosing the mode based on your skin's conditions. Each mode is timed to last for 10 minutes, with a suggested daily total exposure of no more than 40 minutes.

## IV. Care and maintenance

### I . Cleaning

- 1) Before cleaning the device, turn it off and leave it to cool down.
- 2) Wipe the surface of the device with a non-alcoholic wet wipe or soft tissue. Note: If water doesn't do the trick, you may use a cotton swab or a clean cloth dipped in 75% isopropyl alcohol.
- 3) Do not use a steel wire ball or corrosive detergents as they may damage the device.
- 4) Do not clean the device or any accessories in the dishwasher or under running water as they are not waterproof.

### II . Storage

- 1) Keep the device dry. Do not immerse the device or any of its parts in water to avoid the risk of electric shock.
- 2) To prevent damage or loss of the device and its accessories, store them in the packaging box after each use. Before storage, make sure to clean them thoroughly.
- 3) If the device is not used for a long time, store it in a dry, cool place away from sunlight. Avoid contact with hard objects to prevent damage to the device.
- 4) Do not store the device and its accessories in a humid environment for an extended period, as this may affect the performance of the device.
- 5) Regularly check the device (including the cable) for any obvious signs of damage. In case of damage or breakage, stop using the device. To avoid the risk of entanglement or tripping, always unplug the device when not in use.

### III. Storage and transportation requirements

- Ambient temperature: -10-60°C
- Relative humidity: 10%-90%RH
- Atmospheric pressure: 50kPa~106kPa
- Indoors with non-corrosive gases and good ventilation.
- Transportation is subject to the terms stipulated in the purchase contract. However, precautions must be taken to prevent severe impacts, vibrations, as well as exposure to rain or snow during the transportation process.

### IV. Troubleshooting

If you encounter any problems while using the device, you can refer to the frequently asked questions (FAQ) below. If the FAQ does not address your concerns, please contact the manufacturer or local distributor at [seasupport@ulike.com](mailto:seasupport@ulike.com).

## V . FAQ

**Q: Can I use it in the bathroom?**

A: No, This product is not waterproof, so do not use it in the bathroom or near water sources.

**Q: I feel a warm sensation on my face after use. Is that normal?**

A: Yes. The warm sensation is a normal effect of the light therapy. It means that the device is stimulating your skin cells and improving your blood circulation.

**Q: Why does my skin feel dry after use?**

A: Light therapy may cause your skin to become dehydrated. You can use a moisturizing skincare product to hydrate your skin after using the device.

**Q: Can I wash it with water?**

A: No. This product is not waterproof, so do not wash it with water as it may lead to a malfunction of the device.

**Q: How often and how long should I use the product?**

A: You can use the device up to 6 times a week, depending on your skin's condition. Once your skin gets used to the device, you can use it every day. It is recommended to use the device in the morning before applying makeup or in the evening after cleansing. You can choose the mode that suits your needs. The device will automatically shut off after 10 minutes in each mode, and you can use different modes within this period if desired.

※ To prevent excessive care and damage to the skin barrier, you should not use the device for more than 40 minutes a day.

**Q: Can I wear makeup while using the LED mask?**

A: No. You should remove your makeup and cleanse your face before using the LED mask. This will ensure that your skin is free of oil and dirt, and that the light can penetrate your skin more effectively..

**Q: Can it be taken on high-speed trains or airplanes?**

A: Yes. The device is safe and portable, and you can take it with you on high-speed trains or airplanes. However, you should follow the security regulations and instructions of the transportation authorities.

**Q: Does it have any impact on the eyes?**

A: No. The light emitted by the LED mask has no effect on eyesight and is harmless to the human body. Additionally, the eyewear can effectively prevent discomfort caused by light exposure to the eyes.

## VI. Product specifications

Product name	ReGlow LED Face Mask
Product model	UM10
Rated input	5V $\pm$ 1A
Dimensions	235.9*191*145mm
Battery capacity	2600mAh
Charging time	3h
Battery life	30-40min
Mode	Skin brightening mode Skin firming mode Rejuvenation mode Anti-acne mode
Red light wavelength	630nm $\pm$ 10nm
Infrared light wavelength	830nm $\pm$ 10nm
Yellow light wavelength	590nm $\pm$ 10nm
Blue light wavelength	415nm $\pm$ 10nm
Manufacturer	Shenzhen Ulike Smart Electronics Co., Ltd.
Address	810, Building 1, Xunmei Science and Technology Plaza, No. 8 Keyuan Road, Science Park Community, Yuehai Sub-District, Nanshan District, Shenzhen 518000, Guangdong, P.R. China

## VII. Environmental conditions

Storage and transportation conditions		Operating conditions	
Storage temperature	-10~60℃	Operating temperature	10~40℃
Humidity	10%~90%RH	Humidity	15%~80%RH
Relative pressure	50KPa~106KPa	Relative pressure	70KPa~106KPa

## VIII. Warranty terms

- The warranty period for this product may vary by the country of sale, and we will strive to extend it where permitted by local laws.
- The warranty period for this product is 1 year from the date of purchase.
- In the event of any quality issues within the warranty period, Ulike will repair or replace the defective parts at no cost. However, damage caused by improper use or non-compliant operations is not covered. In such cases, Ulike gives you the option to either have damaged parts replaced at the original price or opt for a trade-in. The ownership of the replaced parts or old device shall transfer to Ulike.
- Upon request for service, the warranty certificate must be presented. The warranty period for Ulike products starts from the date of purchase. If proof of purchase cannot be provided, the manufacture date as recorded by Ulike will be deemed to be the start of the warranty period.
- Ulike shall not reissue purchase invoices or warranty certificates. Please keep them properly.
- In-warranty repair or replacement may not apply under the following circumstances:
  - 1) Damage caused by unpredictable events such as fire, earthquakes, or floods.
  - 2) Failures or damage resulting from your negligence, misuse, or use under abnormal conditions.
  - 3) Damage caused by disassembly without authorization or improper storage (such as damage by rodents or liquid ingress).
  - 4) Damage caused by the use of the product in non-standard environmental conditions, including but not limited to the use of the product beyond the recommended workload.
  - 5) Damage caused by situations beyond Ulike's control.
- Ulike does not assume responsibility for consequential damage resulting from the functionality or use of this product, including but not limited to loss of business profits, business interruption, other losses caused by the use of this product, costs incurred in replacing equipment, devices, or services due to the use of the product, and any claims from your customers.
- For the type and extent of liability arising from unforeseen events, Ulike will adhere to applicable laws and regulations.

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