
# Ulike

## Air 10 IPL Hair Removal Device User Manual

Air Series

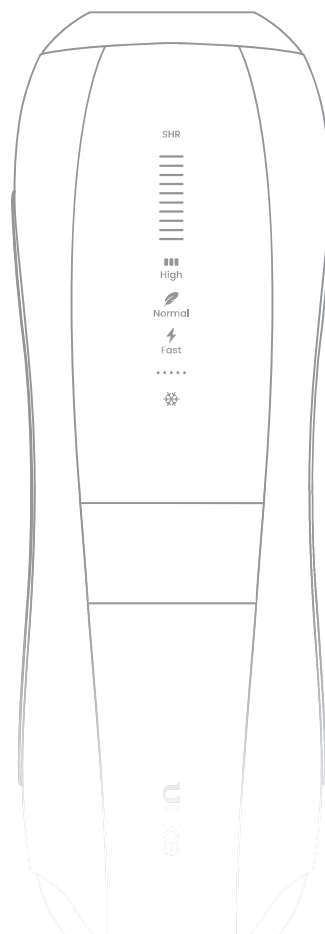ENGLISH ————— 02

Hebrew ————— 22

Thank you for choosing Ulike.

Founded in 2013, Ulike has focused on at-home IPL hair removal devices from day one.

We believe hair removal should be a safe, smooth, and pain-free experience.

And this device is going to be the hair removal system you like.

Feel free to contact us via email, telephone, or chat with us through our website if you need assistance or have any questions. We're here to help!

EU: [de.ulike.com](http://de.ulike.com) & [de.support@ulikeglobal.com](mailto:de.support@ulikeglobal.com)

UK: [uk.ulike.com](http://uk.ulike.com) & [uk.support@ulikeglobal.com](mailto:uk.support@ulikeglobal.com)

FR: [fr.ulike.com](http://fr.ulike.com) & [jpsupport@ulikeglobal.com](mailto:jpsupport@ulikeglobal.com)

EU Amazon: [eusupport@ulike.com](mailto:eusupport@ulike.com)

Shenzhen Ulike Smart Electronics Co., Ltd.

810, Building 1, Xunmei Science and Technology Plaza, No. 8 Keyuan Road, Science Park Community, Yuehai Sub-District, Nanshan District, Shenzhen 518000, Guangdong, P.R. China

Before using the device, please read the user manual carefully, including all contraindications, warnings and safety information, and keep them for future reference. Illustrations in this manual are for reference only.

## Safety Precautions

### ■ Contraindications

Not everyone is suitable to use the IPL Hair Removal Device, please refer to the following information in detail before using the device.

- DO NOT use the device on skin including burns, tattoos, or scars to avoid injury.
- Do NOT use the device if your skin has been artificially or naturally tanned recently or burned from overexposure to sunlight. Your skin may be extra sensitive following exposure and particularly susceptible to side effects of IPL treatment (e.g. burn, blister, discolor, or scar your skin).
- DO NOT use if you are sensitive to sunlight (photosensitivity) or are taking medication that causes photo sensitivity e.g., certain antibiotics such as tetracycline, NSAIDS such as ibuprofen or retinoids such as Accutane and/or topical retinoids such as Retin A. Always check the instructions that come with your medical products to see if photosensitivity is listed as a side effect. Seek for medical advice if you are unsure of whether it is suitable to use IPL Hair Removal Device. Using the device on photosensitive skin could result in skin damage including blisters and swelling.
- DO NOT use the device if you have a history of skin cancer or pre-cancerous lesions (e.g. nevi or a large number of moles).
- DO NOT use if you are pregnant, lactating or menstruating.
- DO NOT use the device on people under the age of 18. The effects of using the device in this age group are unknown.

### ■ Warnings

- DO NOT use the device on the following areas:
  - On testicles, genitals and nipples.
  - On dark brown or black spots such as birthmarks, moles or warts.
  - On tattoos or long-lasting make up.
  - Do not use the device on areas affected by skin conditions such as psoriasis, vitiligo, eczema, acne, herpes simplex, active infections, or wounds.
  - On areas that have had plastic surgery.
  - On areas near eyes.
  - Men should not use this product on the testicles and penis.
- DO NOT exceed the usage specified in this manual. Damage of the skin may occur after prolonged or repeated surface application on one site.

### EYE SAFETY

The device emits intense pulses of light during treatment, and direct exposure can cause damage to the eyes.

- DO NOT look directly into the treatment window when the device has electrical power and is switched on.
- DO NOT attempt to activate the device towards the eyes.
- DO NOT treat areas around the eyes (eyebrows or eyelashes).

It is recommended to wear goggles and look away from the device when activating the lighting button. Additionally, it is advised to avoid having the air inlet facing toward your eyes during treatment. This allows for proper contact of the device with the treatment area, ensuring effective light transmission while shielding the eyes from stray light.

### ELECTRICAL AND FIRE SAFETY

The PATIENT is an intended OPERATOR.

- DO NOT disassemble, repair and refit the product without presence of a professional personnel as this may cause damage, malfunction or personal injury.

After-service is provided by an authorized service spot. If needed, circuit diagrams, component part lists, descriptions, calibration instructions can be provided by MANUFACTURER to assist to SERVICE PERSONNEL in parts repair.

Any unauthorized modification is prohibited, which may lead to unacceptable risk,

- DO NOT use if the device is damaged or broken; for example, cracked housing, cable damage (internal wires visible), cracked or broken treatment window etc.
- DO NOT use near bathtubs, showers, basins or other vessels containing water. The device should not be used if it is exposed to water.
- DO NOT use if the device becomes uncomfortably hot to touch.
- DO NOT attempt to activate (flash) the device onto any surface other than skin.
- DO NOT block the air vents on both sides of the device.
- DO NOT carry out maintenance during use.
- DO NOT position the IPL Hair Removal Device in a location away from the electrical outlet.
- DO NOT leave the device unattended when it is switched on to avoid the risk of fire or burns.

### ■ Unsuitable Situation

- This device is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device to avoid the risk of fire or burns.
  - This device is not intended for use by children. Children should be supervised to ensure that they do not play with the device to avoid the risk of fire and burns.
  - If the mains cord is damaged, you shall have it replaced by Shenzhen Ulike Smart Electronics Co., Ltd., a service center authorized by Shenzhen Ulike Smart Electronics Co., Ltd. or similarly qualified persons in order to avoid risk of electrocution.
  - Water and electricity are a dangerous combination! To avoid risk of electrocution.
    - do not use this device in wet surroundings (e.g. in the bathroom or near a shower or swimming pool);
    - do not let water run into the appliance.
  - Check if the voltage indicated on the device corresponds to the local mains voltage before you connect the device to avoid risk of electrocution or long-lasting damage to the device.
  - If the wall socket used to power the device has poor connections, the plug of the device becomes hot. Make sure you plug the device into a properly installed wall socket to avoid the risk of fire and burns.
- 
- DO NOT subject the device to heavy shocks to avoid risk of damage to the lamp.
  - CAUTION: Always unplug the appliance after use and in case of a power failure to avoid risk of damage to the device.
  - The risk group classification: Exempt Group.
  - Keep the device away from lint, dust and light, or the performance will be affected.
  - Keep the device out of children and pets' reach.
  - DO NOT expose treated areas to the sun. Wait at least 48 hours after treatment before exposing treated skin to direct sunlight. Your skin may be extra sensitive following IPL treatment and particularly susceptible to sunburn. Apply sunscreen (SPF 15 or greater) to the treated area or cover the treated area with suitable clothing.

### ■ Precautions

- The device is intended for a single user.
- Professional skills and training are not required for the usage, please check the user manual before use.
- This device is not recommended for use on red, grey, blonde or very light hair since it is not effective on these hair colors.

## Product information

Product name: **IPL Hair Removal Device**

Model: UI20 DB, UI20 RE, UI20 GP, UI20 PW

### ■ Indication for Use

IPL Hair Removal Device is indicated for the removal of unwanted hair. The device is also indicated for the long-lasting reduction in hair regrowth, defined as the long-term, stable reduction in the number of hairs regrowing when measured at 6, 9 and 12 months after the completion of a treatment regime.

### ■ Packing list

Upon opening the package, you can see the following contents:

IPL Hair Removal Device\*1

Power Adaptor\*1


Power Cord\*1

User Manual\*1

Razor\*1

Goggles\*1

## ■ Product Structure
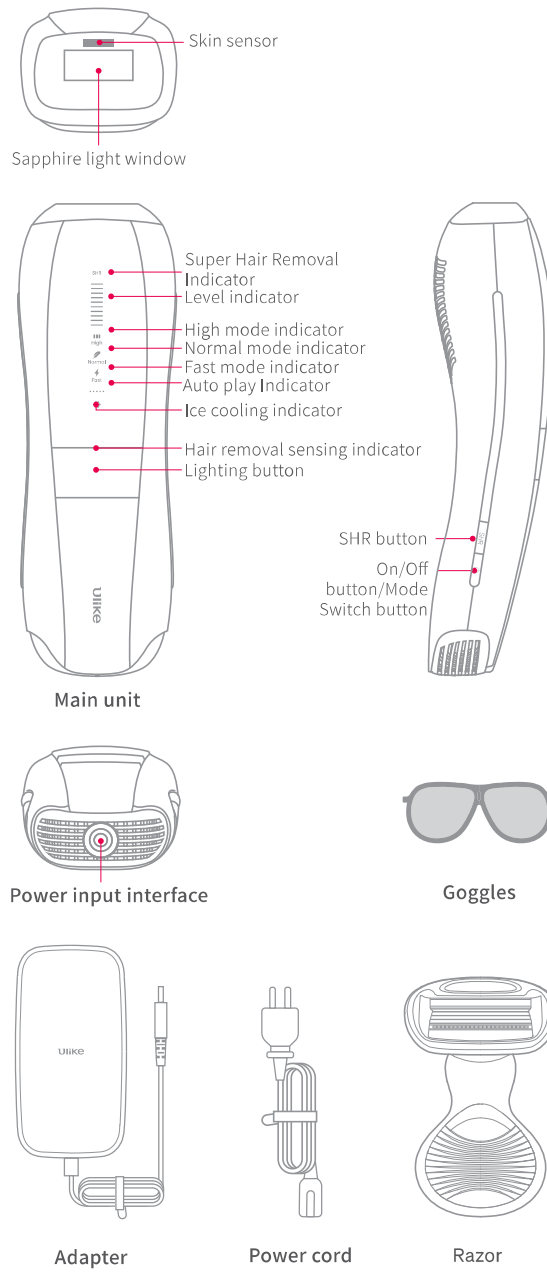
## ■ Operation Principle

IPL Hair Removal Device is a light-based device for long-term hair removal, designed for home use. It is intended for the removal of unwanted hair and long-lasting reduction in hair regrowth. Ideal body areas include the armpit, bikini line, arms and legs. The device uses the IPL technology with lower energy level, including 4 power modes. Intense Pulsed light technology is able to achieve long-term hair removal results at a fraction of the energy level used in other light-based hair removal devices.

The principle of operation is selective photo thermolysis and the mechanism of operation is to disable hair growth using light to preferentially heat the hair bulb.

## How to Use

### ■ Application Ranges

The IPL Hair Removal Device is for hair removal of large areas (e.g. arms, legs, chest) and small areas (e.g. fingers, lip). Please check the Skin Tone Table below to determine whether IPL is suitable for you.

#### USE ON ACCEPTABLE SKIN TYPES AND HAIR TYPES

This device is to only be used on Fitzpatrick Skin Types I-V. Please refer to the chart below.

BEFORE USING THE DEVICE, identify your skin tone on the basis of Fitzpatrick Phototyping Scale, which was developed by Harvard Medical School dermatologist, Thomas Fitzpatrick, M.D., PhD, to classify a person's complexion and their tolerance of sunlight.

- Please IDENTIFY your skin tone and check the following diagram to see whether the product is applicable:

Skin color Hair color	Skin color					
	White	Beige	Light brown	Medium brown	Dark brown	Brown black and darker
White/Grey	✗	✗	✗	✗	✗	✗
Red	✗	✗	✗	✗	✗	✗
Light blond	✗	✗	✗	✗	✗	✗
Light brown	✓	✓	✓	✓	✓	✗
Brown	✓	✓	✓	✓	✓	✗
Dark brown	✓	✓	✓	✓	✓	✗
Black	✓	✓	✓	✓	✓	✗

(Note: X means not applicable, ✓ means applicable)

#### NOTE:

- Choose the comfortable treatment mode according to your skin tolerance and change from fast mode to powerful mode gradually. If there is no adverse reaction, you could change the treatment mode.
- The color of your skin may change, for example because your skin has been tanned: it may then be darker than the skin tone during the

previous treatment. Therefore, it is important that you check the settings to be used before each treatment and set different settings according to different treatment areas. Some areas may be darker than others.

- User with skin type V should use the device with a particular caution because you are more likely to experience adverse events. If you have trouble identifying yourself as Skin Type V versus VI, please consult your physician.


### ■ Prior to Use

Use the razor provided along with the device to remove all visible hair on the skin before using the device. Make sure there is no hair remaining above the surface of the skin as this can cause treatment discomfort. Remove any cosmetics, lotions and cream from the area you wish to treat. Cleanse the treatment area and pat dry.


### ■ Skin Test

- Before your first treatment, we RECOMMEND you test your skin for reactions.
- Select part of body you intend to treat. Connect power adapter and switch on the device. Press sapphire treatment window firmly on the selected skin until the hair removal sensor indicator lights up. Move the device to another test position and change treatment mode, press the Flash Button to do another test treatment. If it feels comfortable, change treatment mode and try it on another test position.
- Do not apply more than one flash on the same spot. Wait 2 hours after testing, to ensure your skin is suitable for treatment and there is no adverse reaction to the light energy. If there is no skin abnormality after 2 hours, and you didn't feel intense hotness during the treatment, then the energy level you used is suitable for your skin. You can continue to test the next energy level till you find the maximum acceptable level for the intended skin area.


### ■ Instructions for Use

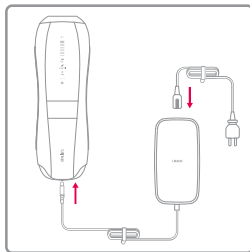
1. Connect the power cord, power adapter and device as shown in figure.

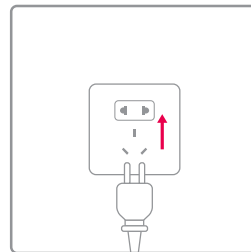2. Plug the other end of power cord into an electrical outlet. Connect the power cord to the power source. (Note: It is a plug-in device, please always keep it connected to the power during the treatment. Otherwise, it will not work).

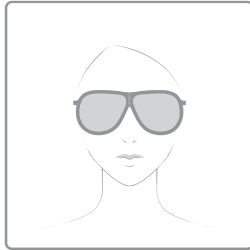
3. Wear the goggles before treatment starts.

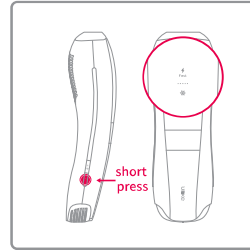4. Short press the power button/mode button on the side of the device to turn it on. The default mode is Fast Mode, and the cooling function will be active during the whole process, so the relevant indicators will be on when the device is turned on.

**Note:** The cooling function is designed for a better hair removal experience, but the function itself doesn't have any therapeutic effect.

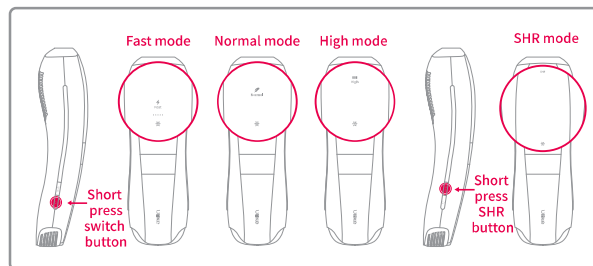5. The device has three modes. When the device is on, short press the power button/mode button to switch to the hair removal mode you need. The mode will change in the following order: Fast Mode - Normal Mode - High Mode - Normal Mode - Fast Mode.


**Note:** Super Hair Removal mode ( SHR Mode )

Area: Armpit and bikini area (coarse hair area)


Recommended Flashes on Each Treatment Area: 2-3 times

**Instructions:**

When starting to use, please try the FAST mode for a skin test first. After confirming that there are no abnormalities in the irradiated area, you can use this product with confidence. The order of use is Fast—Normal—High—SHR. Due to the high energy of SHR, please use the three basic hair removal modes first, and then choose the SHR mode for hair removal after your skin has established tolerance.

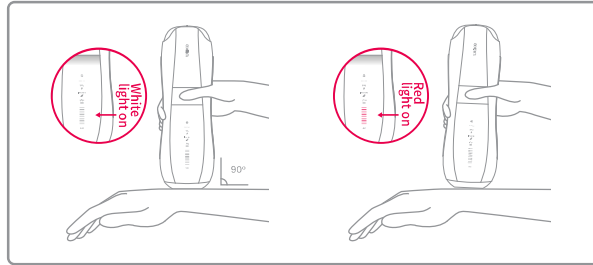6. Place the sapphire treatment window on the skin to be treated vertically and the hair removal sensor indicator will light up in white to indicate that the treatment window is in full contact with the skin. Press the flash button to start hair removal. The red light turning on indicates, which light up in red means that the sapphire treatment window is not in full contact with the skin, and the device will not emit light pulse.

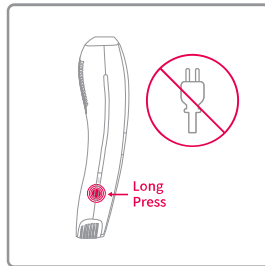
After completing the treatment, press and hold the power button for 1.2 seconds to start the shutdown process. The indicator light turning off indicates that the device is shutting down. Please note that during the shutdown process, the fan will continue to run for 10 seconds to dissipate heat (there will be sound from the fan). The device will only shut down completely after the fan stops spinning (the fan noise disappears). Do not unplug the device directly without shutting it down first.

Note: The device has a fan speed detection module which detects the fan speed. When the fan speed is lower than 70% of the rated speed, the device will automatically shut down, so as to avoid the local high temperature leading to burns.

## ■ Skin Color Recognition

The skin tone sensor can detect and identify the color of the skin, and determine the required intensity based on the recognized skin color. Make sure the skin tone sensor is in full contact with the skin. If a valid skin color is detected, the energy indicator light will turn white; if it is not in full contact with the skin, the energy indicator light will turn red. Under the basic hair removal mode, the intensity of Normal and High modes will depend on the user's skin color. Fast mode only has a continuous mode and only one intensity. Normal and High modes can be selected as single or continuous depending on individual needs.


## ■ Operating Modes

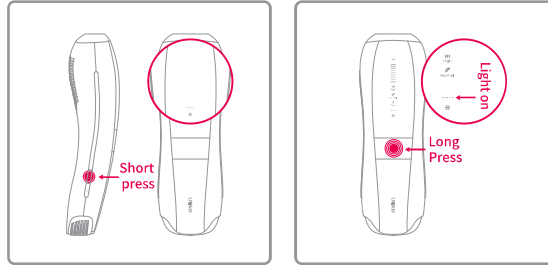The device features two operating modes: auto mode and manual mode. In automatic mode, the device automatically emits light pulses at regular intervals without requiring manual operation. In manual mode, the user must manually press the flash button to emit a light pulse. It's important to note that the SHR mode does not include an auto mode.

### • Auto Mode (for large areas such as arms, legs):


Long press the flash button for 1 or 2 seconds to active the Auto Mode and the auto play indicator will light up. When entering the auto mode, flashes will be emitted at a certain interval, so please keep the sapphire treatment window in full contact with skin and keep moving the device, slide it back and forth twice on the area to be treated to ensure that the area is covered evenly and avoid missing light. Under the auto mode, press the flash button for 1 or 2 seconds again to switch back to the single play mode.

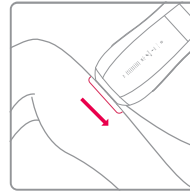**Tip:** You can adjust your movement speed according to the flashing speed, so that the product can cover the area more evenly. The flashing frequency is higher in Fast Mode and slightly lower in High Mode. Furthermore, the higher the level, the slower the light output speed.

### • Manual mode (suitable for small areas such as lip hair):

Place the sapphire treatment window on the area to be treated and fully in contact with the skin. The hair removal sensor indicator will light up, indicating that the treatment window is in full contact with the skin. Press briefly the flash button to start hair removal. Move the sapphire treatment window to the next area, wait for the white light of the hair removal sensor indicator to light up, and then press the flash button again to repeat the above operation in order.

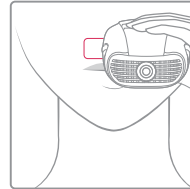**Tip:** If your skin begins to feel excessive heat during the hair removal process, pause for a few seconds to let the sapphire treatment window cool down before using the product again.

### ■ Suggested Treatment Regime

Week 1&4: Three times a week, every other day.

From week 5: Adjust the treatment schedule to once every two weeks or once a month, based on your hair growth.

Please note that you can customize the treatment schedule and level according to your hair type. For thicker or heavier hair, you can consider increasing the level and frequency within your tolerance. Fast Mode is the gentlest, Normal Mode is intermediate, High mode has relatively stronger energy output, and SHR Mode is the most powerful with the highest energy output. Keep in mind that individual hair growth varies, and results may differ for each person.

### ■ Skincare After Treatment

After hair removal, it is recommended to apply body lotion or cold compress gel for skin care. Do not use skin care products containing alcohol. Do not expose treated areas to the sun. Wait at least 48 hours after treatment before exposing treated skin to direct sunlight.

Do not use skincare products containing alcohol or acidic/alkaline substances.

## Maintenance

### ■ Cleaning

In order to ensure the best results and service life of the product, please clean the product before and after each use to avoid debris accumulation.

- Before cleaning, the device must be switched off and unplugged to allow the device to cool.
- Use a soft, clean and dry cloth to wipe and clean the treatment window surface.
- Never wash the product and accessories under faucets or in dishwashers.
- Do not use abrasive cleaners (toothpaste, etc.) or corrosive liquids (gasoline, motor oil, acetone, etc.) to clean the product.

Note: If the water does not clean well, use a cloth dipped in 75% concentration alcohol to clean the device.

When cleaning is finished, please conduct a visual inspection to check whether there is residual dirt. If the device is not visually clean, please repeat the above cleaning steps.

### ■ Storage

- Clean the device before you store it.
- Store the device in a dry, dust-free environment at -10°C ~+55°C / 14 °F ~ 131 °F

## Disposal

In order to minimize hazards to health and the environment and ensure that materials can be recycled, this product should be disposed of at a separate collection facility for waste electrical and electronic device according to local environmental requirements.

## Troubleshooting

Below summaries the most common problems you could encounter with the device. If you are unable to solve the problem with the information below, please contact your local distributor or the seller

Problem	Possible Reason	Solution
Unable to boot normally	Power is not connected	Check if the device is plugged in.
Unable to emit light	Whether the temperature of the treatment window and the body is too high; Whether the treatment window is close to the skin	Wait for the device to cool down before using it; Hold the device close to the skin until the skin contact indicator lights up.
The work indicator is off	Faulty device contact	Check if the power supply is working properly; Repair or return product immediately.
Odor appears when flashing	Filter is contaminated; Hair residue on body surface	Wipe filter; Shave hair and reduce energy level.
Unbearable pain during hair removal	Did not shave off the hair to be removed	Before using this product, please shave off the hair to be irradiated. If shaving causes skin irritation, cut you hair as short as possible and use it in a comfortable position.
	Energy level is too high	Reduce energy level to a comfortable level
	The treatment window is damaged	Repair or return product immediately
After using, the hair removal area turns red	A slight reddening is normal and harmless to the body, and it will soon fade	Use aloe vera gel to soothe the skin or use treatment window for ice soothing.
Unsatisfactory hair removal effect	You use an energy level that is too low for you	Please choose a higher level next time
	Missing flash areas on the hair removal area	Flashes in close proximity and there should be some overlap between the flash areas.
	This device has no effect on your hair color.	If your hair is pale blonde, gray, red or white, the treatment will not be effective.
	Insufficient frequency of use	To successfully remove all hair, we recommend that you follow the recommended treatment regimen.
	Your response to light exposure is slow	This product should be used continuously for at least 5 months, as it can still suppress hair regeneration during this period.
The product becomes hot and humid during use	May be due to long-term device operation	Turn off the device and let it stand, and then let it cool down before using it.

## Technical Specifications

Product Name	IPL Hair Removal Device
Product Model	UI20 DB, UI20 RE, UI20 GP, UI20 PW Note: Model difference only exists in enclosure color, UI20 DB is Purple, UI20 RE is Pink, UI20 GP is Gold, UI20 PW is Silver
Dimensions (W x H x D)	206.73mm*68.68mm*54.29mm
Power Input	24V $\overline{\text{---}}$ 3.75A, 90W
Energy Level	1-10 level
Spot Size	3.9cm <sup>2</sup> /0.6 sq in
Wavelength	550nm-1200nm
Output Energy	5.0J-26.0J
Energy Density	1.28J/cm <sup>2</sup> -6.66J/cm <sup>2</sup>
Pulse Duration $\pm$ 20%	0.4ms-2.0ms
Software Version:	V1.0
Operating conditions	Operating temperature: 5°C ~ 30°C/ 41°F ~ 86 °F Humidity: 5%- 90%RH Relative pressure: 80kPa -106kPa
Storage and transportation conditions	Storage temperature:-10°C ~ 55°C/ 14 °F ~ 131 °F Humidity: 5%- 90%RH Relative pressure: 50kPa-106kPa
Lifetime	300,000 flashes.
Adapter	Manufacturer: ASIAN POWER DEVICES INC. Model: DA-90Q24 INPUT: 100V-240V~, 50/60Hz, 1.5A OUTPUT: 24.0V $\overline{\text{---}}$ 3.75A, 90W
	Manufacturer: Mass Power Electronics Inc. Model: E090-IE240375M2 INPUT: 100V-240V~, 50/60Hz, 1.5A OUTPUT: 24.0V $\overline{\text{---}}$ 3.75A, 90W
After take out the device from the lowest storage ambient temperature (-10°C), the time required for the device to return to normal operation at ambient temperature of 20°C.	2 hours
After take out the device from the highest storage ambient temperature (60°C), the time required for the device to return to normal operation at ambient temperature of 20°C.	2 hours

## EMC & Biocompatibility

### EMC guidance

- **Warning:** Don't use near active HF surgical device and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.
- **Warning:** Use of this device adjacent to or stacked with other device should be avoided because it could result in improper operation. If such use is necessary, this device and the other device should be observed to verify that they are operating normally.
- **Warning:** Use of accessories, transducers, and cables other than those specified or provided by the manufacturer of this device could result in increased electromagnetic emissions or decreased electromagnetic immunity of this device and result in improper operation.
- **Warning:** Portable RF communications device (including peripherals such as antenna cables and external antennas) should be used no closer than 12 in (30 cm) to any part of the device, including cables specified by the manufacturer. Otherwise, degradation of the performance of this device could result.

## Symbols

No.	Symbols	Description
1		Indoor Use
2		"WEEE (Waste Electrical and Electronic device)". The waste products should be handled legally.
3		Keep dry
4		Manufacturer
5		Warning symbol for use with IPL device: Do not use it directly on the eyes, as the intense light emitted by this appliance may cause eye damage.
6		DC pocket
7		Maximum number of identical transport packages/ items which may be stacked on the bottom package, where "5" is the limiting number.
8		Caution

## Warranty

After-sales protection policy fully complies with relevant sales laws and regulations.

From the date of purchase, the product is covered by a 90-Day Return and 2-Year Warranty. However, warranty may be refused under the following circumstances:

- Failure to follow the instructions, operation, maintenance.
- Mistake which causes damage to the product.
- Exceeding the allowable use environment (such as high temperature, corrosion, etc.) causing product damage.
- Non-normal external forces (falling, extrusion, impact, etc.) causing damage to the product.
- Natural force majeure (fire, flood, lightening, impact, etc.) causing damage to the product.
- Use of non-manufacturers accessories and expansion parts causing damage to product.

## Contact information

Shenzhen Ulike Smart Electronics Co., Ltd.  
810, Building 1, Xunmei Science and Technology Plaza, No. 8  
Keyuan Road, Science Park Community, Yuehai Sub-District,  
Nanshan District, Shenzhen 518000, Guangdong, P.R. China

Feel free to contact us via email, telephone, or chat with us through our website if you need assistance or have any questions. We're here to help!

EU: [de.ulike.com](mailto:de.ulike.com) & [de.support@ulikeglobal.com](mailto:de.support@ulikeglobal.com)

UK: [uk.ulike.com](mailto:uk.ulike.com) & [uk.support@ulikeglobal.com](mailto:uk.support@ulikeglobal.com)

FR: [fr.ulike.com](mailto:fr.ulike.com) & [jpsupport@ulikeglobal.com](mailto:jpsupport@ulikeglobal.com)

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