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## Instructions for use


- To be used on clean skin. The best timing for hair removal is right after face wash. It is recommended to use moisturizer or facial lotion before use and to wait for it to dry completely.
- Switch the On/Off to "ON" mode and place the cotton thread tightly to your skin. Move gently against the direction of the growing hair and the hair will be removed. Once finished proceed to the next area. Please be aware that too much pressure will irritate your skin.
- After removing the hair, it is recommended to apply your favorite cooling gel or lotion at removed hair area.
- It is normal if the cotton thread is broken while you are using it.

Step 1: Pass 36cm of cotton thread through the middle hole, keep the wire the same length on both sides.

Step 2: pass the two side threads through the top holes of the steel.

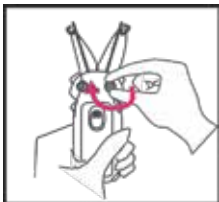Step 3: Entwine the left thread in the bottom of the left bolt and fasten the bolt.

Step 4: Tighten the right thread then entwine in the bottom of the right knob, fasten the bolt.

## Features

Facial threading is a traditional beauty technique from ancient China. By manipulating the threads rapidly, the fine hair on the face can be eliminated. Facial threads will not only leave you with a delicate skin it is also a convenient, safe and fast way to remove unwanted hair. This Hair Remover follows the concept of facial threading to remove unwanted hair as well as constrict pores. The new born hair will be finer as well. Without any chemical ingredients, heat, medicine or knife used, there is no side effect. It is the most ideal hair remover for body and face since it treats your skin in a safe way.

## Warnings


- To avoid malfunction, don't get close to head hair, eyelashes and also avoid contact with metal strings.
- Never leave your hair remover unattended when in use, household use only.
- Use the hair remover only for purposes described in this instruction manual. Only use assorted adapter and supplied cotton thread.
- To avoid malfunction, keep the hair remover and adapter away from water.
- Patients under treatment for varicosity, diabetes, hemophilia or other disease, and people with freckles, Mongolian spots, sensitive skin, and low immunity system including pregnant women should consult your doctor before use.
- If any irritation occurs, stop using the device and consult your doctor.
- Please keep out of children reach.
- For proper result, please keep the hair remover clean.
- It is recommended before using on facial hair or thick hair to first test it on arm hair. Only after this trial you may start the use.
- Do not change the chargeable batteries. It is forbidden and dangerous to use this product with AA or AAA batteries.

Thank you for using our epilator. Please read this manual indetail to ensure better use effect for a more comfortable use. Please keep this manual for future reference.  
If you have problems during use, please contact our customer service: 09-9725151.

## Technical Parameters

Name	Hair Remover	Rated Power	5W±0.1
Charging Voltage	5V±0.2	Charging Time	6-8 Hours
Batteries	2 pcs of 500mAh rechargeable batteries		

## Parts Description
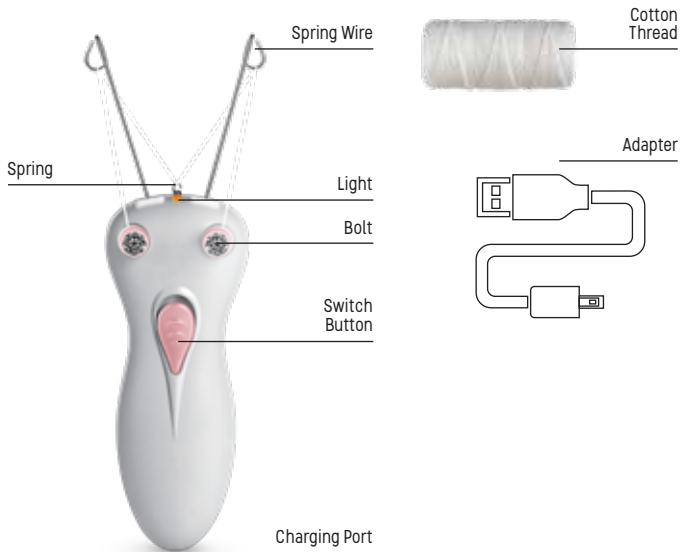
# HAIR REMOVAL THREAD DEVICE

Instructions Manual

Model: SG-667